

A Dedicated Team Just For You



Every Bob and Brad product is backed by a 1-year warranty and an award-winning customer service team. If you are in any way dissatisfied with a product you ordered, we'll exchange it, replace it or refund your money under the warranty. Simply email us the problem, and we'll take care of it.



Troubleshoot

Help you solve any issues with your device.



Product Knowledge

Answer questions about your device and best practices.



Initiate Returns

Provide instructions for replacements, returns and refunds.



Place New Orders

Help with placing orders online and/or authorized resellers.

Get in Touch


Call: +1 (612) 567-3035 Mon-Fri 9:00AM-4:30PM (CST)

Email: support@bobandbrad.com

FLIGE Products LLC

3001 Esperanza Crossing Austin, TX 78758

**Bob
AND
Brad®**

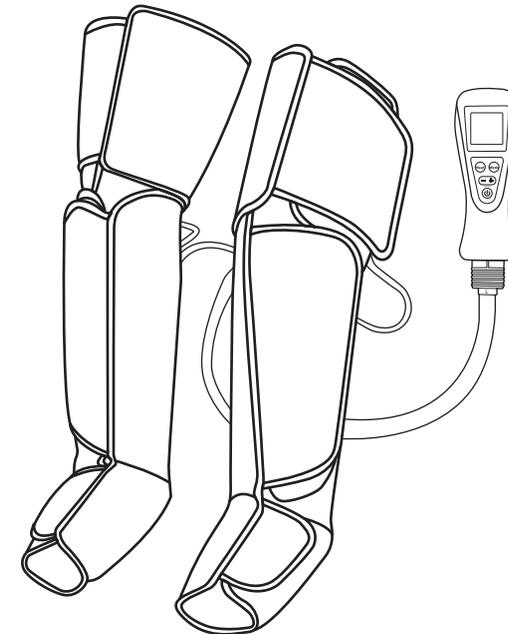

AIR COMPRESSION LEG MASSAGER

User Manual

Model: HY-1146

Product Specifications

Product Name	Air Compression Leg Massager
Product Model	HY-1146
Input	12V == 2A
Power	24W
Operating time	15 minutes

Packing List

1 x Air Compression Leg Massager
1 x Power Adapter
1 x User Manual

FCC Statement

This device complies with Part-15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

FCC Statement NOTE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installation. This equipment generates, uses, and can radiate radio frequency energy, if not installed and used in accordance with the instructions. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause any harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Plug the equipment into the outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help with the supplied fastener. **DO NOT** hang the unit by the power cord.

Danger

1. Do not use spare parts that are not approved by the manufacturer.
2. Non-professional repair personnel should not dismantle and/or replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by the related professional personnel in order to avoid danger.
4. Do not use or store the product in a place with humidity or high temperature.
5. Unplug the power cord after each use.
6. Do not use the product longer than 15 minutes at a time.
7. If any adverse reactions occur, stop use immediately.
8. If the temperature is too high, or the product begins to smell or emit an odor, please stop using it immediately and cut off the power.
9. Do not move the product by pulling the power cord, to prevent the power cord from getting loose.
10. This product is not applicable to those who are deficient in physical, sensory, or mental abilities; or, lack experience and knowledge in this regard (including children).
11. People who are not sensitive to heat should use the heating function with caution.

Cautions


1. If you fall into any of the following categories, please consult a doctor before use:
 - Pregnant women, children, or elderly
 - People who suffer from skin infections or have redness and burns on the skin
 - People who use a pacemaker or similar implantable medical device

- Patients with heart disease or malignant tumors
 - People who are receiving medical treatment or feel discomfort
 - People suffering from osteoporosis
 - People with high blood pressure or abnormal blood pressure
2. This product is not used for any medical purposes. It is designed for personal and non-professional use only.


Cleaning and Maintenance

1. Before cleaning, please turn off the device and unplug the power cord.
2. Wipe the product's surface with a soft, dry cloth.
3. Do not clean the product or its accessories with corrosive detergents.
4. Please store in a cool, dry place when not in use.
5. Avoid exposure to direct sunlight and humid environments. These may cause color fading and/or damage to the product.
6. If the power adapter is damaged, it must be replaced by FLIGE or an authorized manufacturer.

Product Layout

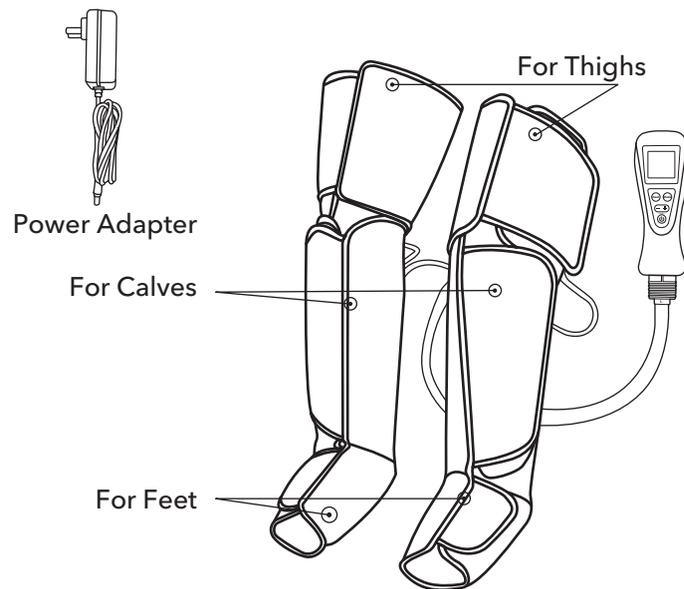

Control Panel

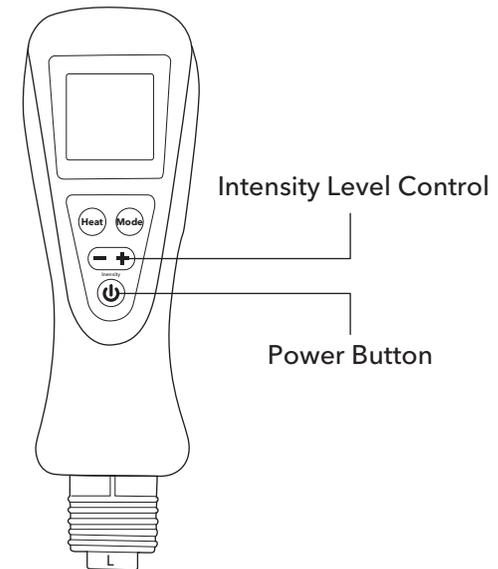

How To Use

Connect the adapter to the power supply and the DC power cord; the product will then be in standby mode.

- ⏻ Press this button when the device is in standby mode; the device will turn on.

Check that all function buttons and their corresponding indicator lights are working properly.

The massager will be in its first massage setting and second intensity level, by default. Once powered on and working, the automatic shut-off countdown timer will begin (15 minutes).

If the button is pressed when the device is in use, or 15 minutes have elapsed, all massaging actions will stop and the system will automatically turn off and return to standby mode.

- 🔥 Press this button to start or switch the heat level.

 Press this button to switch massage areas.

Mode 1: Foot Massage

Mode 2: Calves Massage

Mode 3: Thigh Massage

 Press this button when the device is in its air compression massage setting to increase the intensity. Every time the button is pressed, air compression intensity switches to a higher level (up to level 4). The corresponding indicator light indicates the exact level of intensity.

 Press this button when the device is in its air compression massage setting to lower the intensity back to level 1. The corresponding indicator light indicates the exact level of intensity.
If any button is pressed during the air compression massage setting, the airbags will deflate for 2 seconds. Then, the next massage action will begin.

Troubleshooting

PROBLEM	CAUSE	SOLUTION
The product does not work	1.The power is not connected or the plug is loose or falls out, resulting in poor connection 2. The power switch is not turned on or timer ends	1.Plug in the power cord correctly 2.Turn on the power switch
Some part can't inflate	The current massage mode excludes the area	Please choose the corresponding mode
	The tubes have not been attached securely, and there are air leaks	Please check that the tubes are attached securely to the boots (no gaps between the two)

PROBLEM	CAUSE	SOLUTION
Emit weak sounds when in use	Sounds emitted from the motor during operation.	These are normal operating sounds from the motor
Massage feels too strong or uncomfortable	Legs wrapped too tightly	1.Choose the lowest intensity 2.Open the Velcro of the corresponding part 3.Press POWER button to turn off
	Set the highest intensity	4.Disconnect the corresponding air hose, the airbag of the corresponding part will quickly deflate
The product suddenly stops working	The massager's overheating protection feature will switch off the massager automatically when it's too hot	Wait for the massager to cool down before powering it back on
	15-minute timer finished	Press the power button again
The handheld controller gets hot	The product has been used continuously for too long	Use for max 15-20 minutes, then turn off and allow to cool before next use to avoid OVER USE
Feel a difference of intensity between boots when massaging both legs	A tube has not been attached correctly, or is attached to the wrong hose	Please make sure the tubes are attached securely
	When wearing the boots, the gap you leave may be different	Adjust to leave the same sized gap in both boots
If you encounter any other problems, contact: support@bobandbrad.com		